



BLARH1-10

Amplifier 200-600MHz Operating & Service Manual

Version 001

The information in this manual may be altered without notice.

BRUKER accepts no responsibility for actions taken as a result of use of this manual. BRUKER accepts no liability for any mistakes contained in the manual, leading to coincidental damage, whether during installation or operation of the instrument. Unauthorised reproduction of manual contents, without written permission from the publishers, or translation into an other language, either the entire manual or a part of it, is forbidden.

This manual describes the units as they are at the date of printing. On request, the manufacturer shall supply circuit diagrams, lists of components, descriptions, calibrating instructions and any other information for use by qualified personnel of the user, in charge of repairing the parts of the unit which have been stated by the manufacturer to be "repairable". Such supply shall in no event constitute permission to modify or repair the units or approval of the same.

All rights reserved for the units, circuits, processes and appellations mentioned herein.

This unit is not designed for any type of use which is not specifically described in this manual. Such use may be hazardous.

This manual was written by

[BARTHELEMY Philippe](#)

© April 26, 2005: Bruker Biospin SA

Wissembourg, France

P/N: Z31739
DWG-Nr: 1446.001

Contents

	Contents	iii
1	General information	5
1.1	Introduction	5
1.2	Initial turn and procedure	5
2	Safety	7
2.1	Labels	7
	Identifying plate	7
	Manufacturer's nameplate	8
2.2	Safety labels and symbols	9
	Dangerous area labels	9
3	Installation	11
3.1	Initial Inspection	11
	Mechanical Check	11
	Claim for Damage	11
	Reshipment and Repackaging Requirements	11
	Auxiliary Kit	11
3.2	Installation Requirements	12
	Bench Operation	12
	Cooling and Ventilation	12
3.3	Power Requirements	12
3.4	System Check	12
3.5	Initial Turn On Procedure	12
4	Operation	13
4.1	Front Panel description	13
4.2	Indicators	14
4.3	Connectors	14
	True tables	15
4.4	TCU T1 Signal Accessory Connectors	15
4.5	Rear Panel	17
4.6	BLARH1-10 Amplifier 200-600MHz bloc diagrams	18
5	Technical description	21
5.1	System overview	21
5.2	Theory of operation	21
	RF Amplifier path	21
	RF Outputs router without QNP Relay	22
	RF Outputs router with QNP Relay	23
	Control Board	24

Contents

	Status Board	24
6	Specifications	25
6.1	General specifications	25
	Figures	27
	Tables	29

General information

1

Introduction

1.1

The BLARH1-10 Amplifier 200-600MHz (P/N:W1345023) is a double ^1H band low power 1W & 10W amplifier including RF Output DPDT Switches, specifically suited for use with an external 100W-150W ^1H high power amplifier (such a BLAH100 Amplifier 200-600MHz (P/N: W1345060) or a H150 output from a BLAH1000 Amplifier) to obtain a BLARH100 RF topology.

This amplifier is fitted for ^{19}F decoupling during ^1H observe or ^1H decoupling during ^{19}F observe.


It includes "TCU-T1 Signal" accessory specially suited for the use with an old "solid equipment" (MSL or DSX) to dispatch real time "NMR5_13 & NMR5_14" and "relays" signals.

This accessory allows to dispatch the power of both channels (^1H amplifier channel & X amplifier channel) through three outputs: ^1H ; ^{19}F and, XQNP in option.

Except main power supply, all sockets are coming through the front panel.

Two "BNC" sockets in front panel receive "real time" switching signals.

Figure 1.1. BLARH1-10 Amplifier 200-600MHz without XQNP accessory

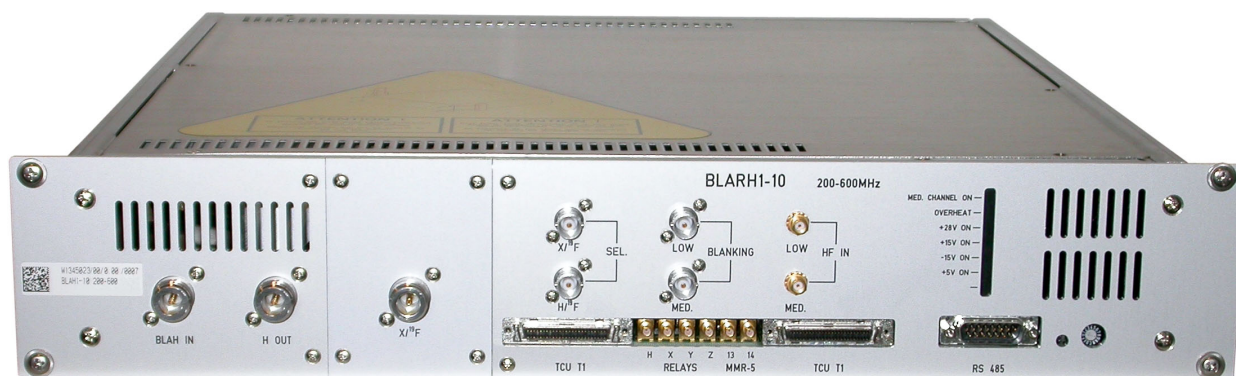

Initial turn and procedure

1.2

Before inserting main power cord, be sure that all RF (HF in and HF out) channels, select signals (SEL. X/ ^{19}F & SEL. $^1\text{H}/^{19}\text{F}$), Blanking and RS485 are correctly plugged.

General information

Labels

2.1


Labels are provided to alert operating and service personnel to conditions that may cause personal injury or damage to the equipment from misuse or abuse. Please read the labels and understand their meaning.

Identifying plate

2.1.1

The BLARH1-10 can be identified by an identifying plate at the front panel of the unit that has following information.


Figure 2.1. Identifying plate

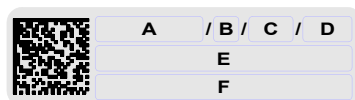

- **(A) Part Number**
This field indicates the part number of the product.
- **(B) Variant**
This field indicates the variant number that identifies the production category of the product. The default variant is 00.
- **(C) ECL**
This field indicates the revision number that identifies the product configuration. The initial revision is 0.00.
- **(D) Serial Number**
This field indicates the serial number of the product.
- **(E) Type**
This field contains the designation of the product.
- **(F) Information**
This field contains additional information about the product.

The BLARH1-10 can be identified by a manufacturer's nameplate at the back panel of the unit that has following information:

Figure 2.2. Manufacturer's nameplate

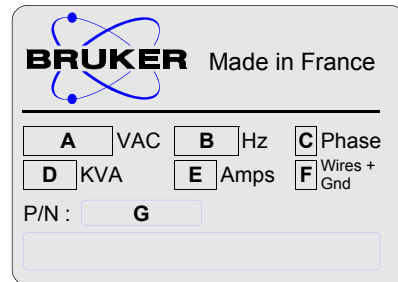

- **(A) Voltage**
This field indicates the input mains voltage of the product.
- **(B) Frequency**
This field indicates the input mains frequency of the product.
- **(C) Phases**
This field indicates the number of phases of the mains.
- **(D) Power**
This field indicates the absorbed power of the product.
- **(E) Current**
This field indicates the absorbed current of the product.
- **(F) Wires**
This field indicates number of wires with the ground in the mains cord.
- **(G) Part Number**
This field indicates the assembly number that identifies the part number of the product.

Safety labels and symbols

2.2

Dangerous area labels

2.2.1



WARNING! Risk of electrical shocks when the device operates.
This device is low voltage supplied (+28V), but high power RF transiting signals generate high voltage.

Figure 2.3. Electrical hazard symbol



Figure 2.4. General hazard symbol



Installation

3

Initial Inspection

3.1

Mechanical Check

3.1.1

If damage of the shipping carton is evident, request the carrier's agent be present when the instrument is unpacked. Check the equipment for damage and inspect the cabinet and panel surfaces for dents and scratches.

Claim for Damage

3.1.2

If the unit is mechanically damaged or fails to meet specifications upon receipt, notify BRUKER or our representative immediately. Retain the shipping carton and packing material for the carriers inspection as well as for subsequent use in returning the unit if necessary.

Reshipment and Repackaging Requirements

3.1.3

Whenever possible, the original carton and packing material should be used for reshipment. If the original packing material is not available, wrap the instrument in heavy paper or plastic. Use a strong shipping container. If a cardboard is used, it should be at least 200 lbs. test material.

Use shock absorbing material around all sides of the instrument to provide a firm cushion and to prevent movement inside the container wall on each side. Protect the front panel by means of cardboard spacers inserted between the front panel and the shipping carton. Make sure that the instrument cannot move in the container during shipping. Seal the carton with a good grade of shipping tape and mark the container:

FRAGILE ELECTRONIC INSTRUMENT

Auxiliary Kit

3.1.4

The BLARH1-10 Amplifier 200-600MHz is shipped with an accessories kit containing following items:

- Manual,
- Accessory, supply cord and RF wire cords set.

Installation

Installation Requirements

3.2

No special precautions are necessary. Mount the equipment in an area which is relatively free of vibration, and has sufficient room for cable connections.

Bench Operation

3.2.1

The units can be placed onto a secure flat surface.

Cooling and Ventilation

3.2.2

No specific cooling or ventilation is required. It should, however, be in an environment which conforms the 0° - 50°C (32°F - 158°F) specification, and in an area that does not obstruct the free flow into and out of the unit.

Power Requirements

3.3

BLARH1-10 Amplifier 200-600MHz is powered by the regular 230VAC main. It taps to the maximum 400VA when it is switched on. The connection to this power supply is realized via supply Cord.

System Check

3.4

Before applying RF power for the first time the following items should be checked:

- The RF output an RF input of the accessory must be wired.
- Two external "real time" signals must drive the both selection inputs (SEL. H/¹⁹F & SEL. X/¹⁹F) only one SEL. H/¹⁹F if the "QNP accessory" is not present.
- Don't omit to wire the two RF outputs from the amplifier to the RF inputs of the accessory (use the two "N" ↔ "N" auxiliary kit cords).

Initial Turn On Procedure

3.5

The following list describes how to turn on the BLARH1-10 Amplifier 200-600MHz and what should be seen as this occurs.

Before starting this procedure, make sure that you have properly followed instructions in the **"System Check" on page 12.**

1. Connect the amplifier to the power supply and turn the circuit breaker, to ON.
2. Observe the indicators on the front panel of the amplifier:
 - The +28V ON LED's will illuminate,
 - The +15V, -15V and +5V ON LED's will illuminate.

System is now fully operational.

Operation

4

Front Panel description

4.1

The BLARH1-10 Amplifier 200-600MHz front panel is provided with:

- "BNC" selection and Blanking connectors,
- "N" RF power sockets and "SMA" RF input sockets,
- "Sub D15" socket for the RS485 communication,
- 6 "SMB" and 2 "SCSI" extra connectors used for TCU-T1 signal accessory,
- Status display.

Figure 4.1. BLARH1-10 Front Panel Design without QNP option

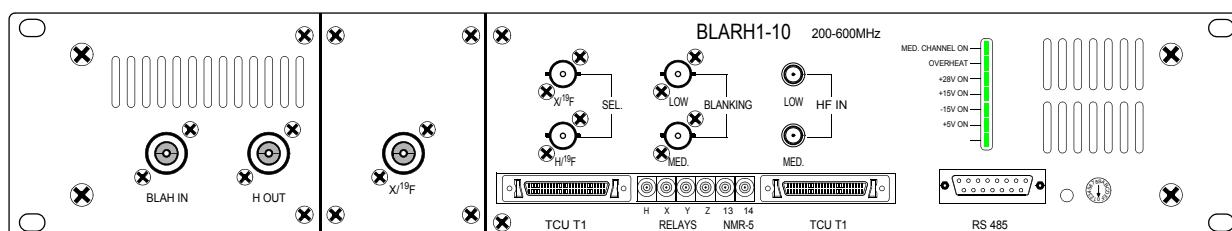


Figure 4.2. BLARH1-10 Front Panel Design with QNP option

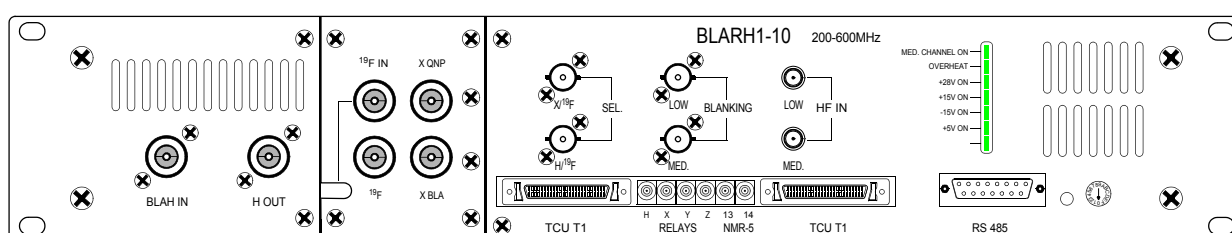


Table 4.1. LED's indicators assignment

MED CHANNEL ON	Indicates that "Med channel" (10W) is sending RF
OVERHEAT	Indicates an overheating or a fan default
+28V ON	Indicates that +28V power supply is ready
+15V ON	Indicates that +15V power supply is ready
-15V ON	Indicates that -15V power supply is ready
+5V ON	Indicates that +5V power supply is ready

Table 4.2. RF Connectors assignment without QNP Relay Option

LOW HF IN	RF in "SMA" Connector (female) 1W channel IN : nominal +4dBm
MED HF IN	RF in "SMA" Connector (female) 10W channel IN : nominal +4dBm
BLAH IN	RF in "N" Connector 100W ¹ H input (for example : BLAH1000 200-400MHz)
H OUT	RF out "N" Connector ¹ H out from Hin + 1W/10W amplifier via DPDT Pin Switch
X / ¹⁹ F	RF out "N" Connector ¹⁹ F out from Hin + 1W/10W amplifier via DPDT Pin Switch

Table 4.3. RF Connectors assignment with QNP Relay Option

LOW HF IN	RF in "SMA" Connector (female) 1W channel IN : nominal +4dBm
MED HF IN	RF in "SMA" Connector (female) 10W channel IN : nominal +4dBm
BLAH IN	RF in "N" Connector 100W ¹ H input (for example : BLAH1000 200-400MHz)
X BLA	RF in "N" Connector 300/500W X input coming from BLAX300 or BLAX500
¹⁹ F IN	RF in "N" Connector X/ ¹⁹ F out from Hin + 1W/10W amplifier via DPDT Pin Switch
H OUT	RF out "N" Connector ¹ H out from Hin + 1W/10W amplifier via DPDT Pin Switch
¹⁹ F	RF out "N" Connector ¹⁹ F out from Hin + 1W/10W amplifier via DPDT Pin Switch and DPDT QNP Relay
X QNP	RF out "N" Connector 300/500W X input via DPDT QNP Relay

Table 4.4. Blanking Connectors assignment

LOW BLANKING	TTL level in "BNC" Connector : 1 = blanked ; 0 = channel active
MED BLANKING	TTL level in "BNC" Connector : 1 = blanked ; 0 = channel active

Table 4.5. RF Output Selection

X / ¹⁹ F	TTL level in "BNC" Connector : 1 = through ; 0 = crossed ^a
H / ¹⁹ F	TTL level in "BNC" Connector : 1 = through ; 0 = crossed

a. For DPDT QNP Relay only

True tables

4.3.1

DPDT Pin Switch ¹H & ¹⁹F outputs routes as a fonction of SEL H/¹⁹F control signal.

Table 4.6. DPDT Pin Switch routing true table

SEL H/ ¹⁹ F	1W (LOW) + external Hin	10W (MED)
TTL Logic 1 (5V)	⇒ ¹ H	⇒ ¹⁹ F
TTL Logic 0 (0V)	⇒ ¹⁹ F	⇒ ¹ H

DPDT QNP Relay BLAH IN & X BLA inputs to ¹⁹F & X QNP outputs routes as a fonction of SEL X/¹⁹F control signal. For more details, see the chapter **"Technical description" on page 21**

Table 4.7. DPDT QNP Relay routing true table

SEL X/ ¹⁹ F	Hin (if SEL H/ ¹⁹ F is LOW)	Xin
TTL Logic 1 (5V) in direct wire mode	supplied to ¹⁹ F	supplied to X QNP
TTL Logic 0 (0V) in reverse wire mode	supplied to X QNP	supplied to ¹⁹ F

TCU T1 Signal Accessory Connectors

4.4

These connectors are defined to supply auxiliaries signals used with old "solid equipment" (MSL or DSX) : please, for more details refer to the Bruker Advanced Service Handbook (BASH : Dual bay > AQX > TCU T1) .

Operation

Table 4.8. TCU T1 Connectors assignment

Pin	Signal	Meaning	Pin	Signal	Meaning
1	word3_0	BLNK_TR1	26	GND	
2	NMR5_0	GAIN0_TR1	27	NMR5_1	GAIN1_TR1
3	NMR5_2	C/AB_TR1	28	GND	
4	word3_1	BLNK_TR2	29	GND	
5	NMR5_3	GAIN0_TR2	30	NMR5_4	GAIN1_TR2
6	NMR5_6	C/AB_TR2	31	GND	
7	word3_4	BLNK_TR5	32	GND	
8	NMR5_6	C/AB_TR5	33	GND	
9	word3_5	BLNK_TR6	34	GND	
10	NMR5_7	RELAY_H	35	NMR5_8	RELAY_X
11	TRIG3		36	GND	
12	NMR5_9	RELAY_Y	37	NMR5_10	RACK_ON/OFF
13	NMR5_11		38	GND	
14	NMR5_12	RELAY_Z	39	GND	
15	NMR5_13		40	GND	
16	NMR5_14		41	GND	
17	NMR5_15		42	GND	
18	NMR6_0	STP1_DIR	43	NMR6_1	LB_SEL
19	NMR6_2	DCM_STRT	44	GND	
20	NMR6_3	STP1_CLK	45	GND	
21	NMR6_4	SPTP2_CLK	46	GND	
22	NMR6_5	RES_STP1	47	GND	
23	NMR6_6	DCM_RES	48	GND	
24	NMR6_7	GO_POS	49	GND	
25	TRIG2		50	GND	

Table 4.9. Correspondance between TCU T1 and SMB Connectors

TCU T1 Connectors	SMB Connectors
NMR5_7 (Pin 10)	Relay H
NMR5_8 (Pin 35)	Relay X
NMR5_9 (Pin 12)	Relay Y
NMR5_12 (Pin 14)	Relay Z
NMR5_13 (Pin 15)	NMR5_13
NMR5_14 (Pin 16)	NMR5_14

Rear Panel**4.5**

The rear Panel of the BLARH1-10 Amplifier 200-600MHz is free of elements in exception of the three pole (2P + E) line socket.

- No power switch required. Power cord insertion will switch on the device.
- No fuse required. There are located in the power supply unit.

Figure 4.3. Amplifier block diagram without QNP option


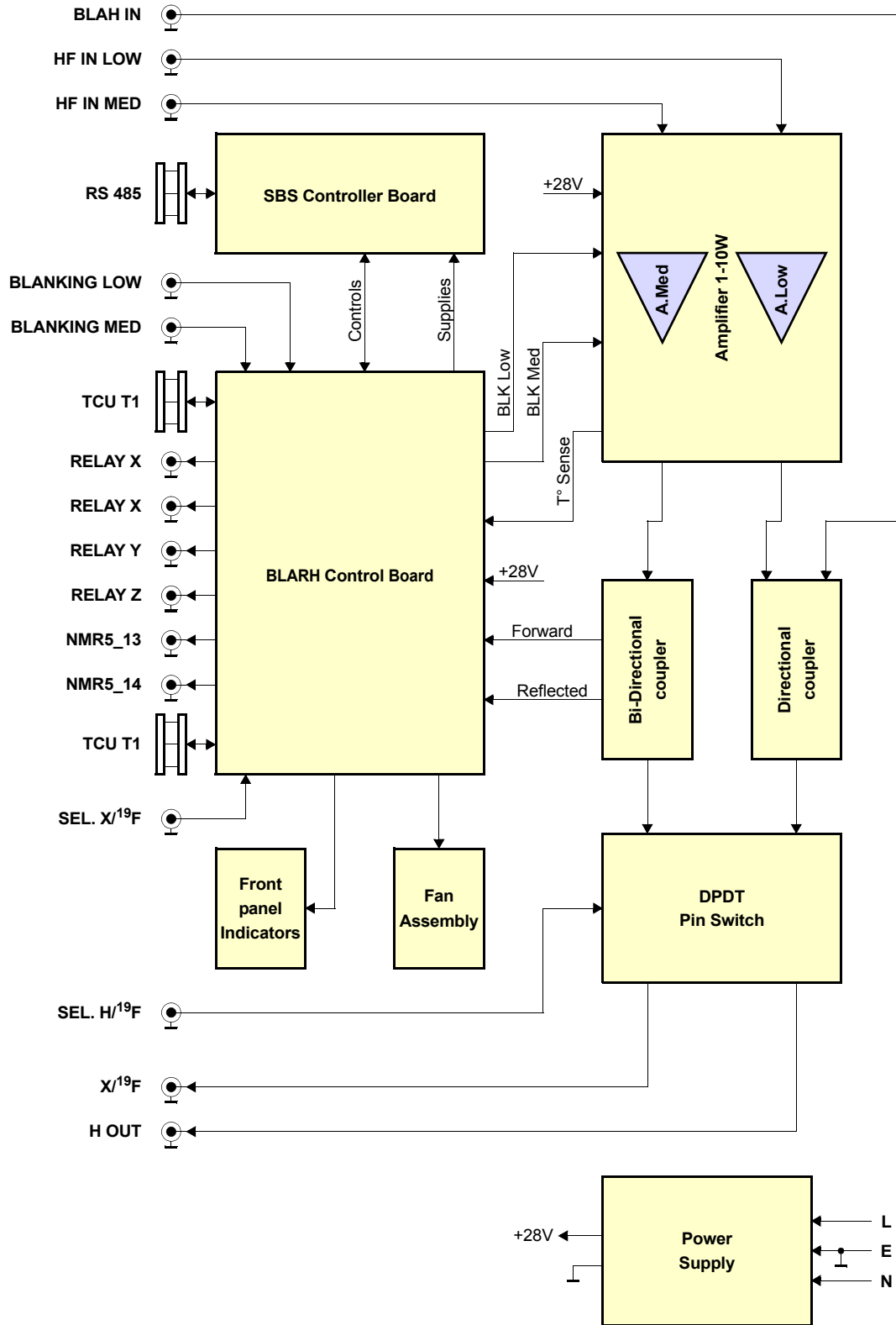


Figure 4.4. Amplifier block diagram with QNP option

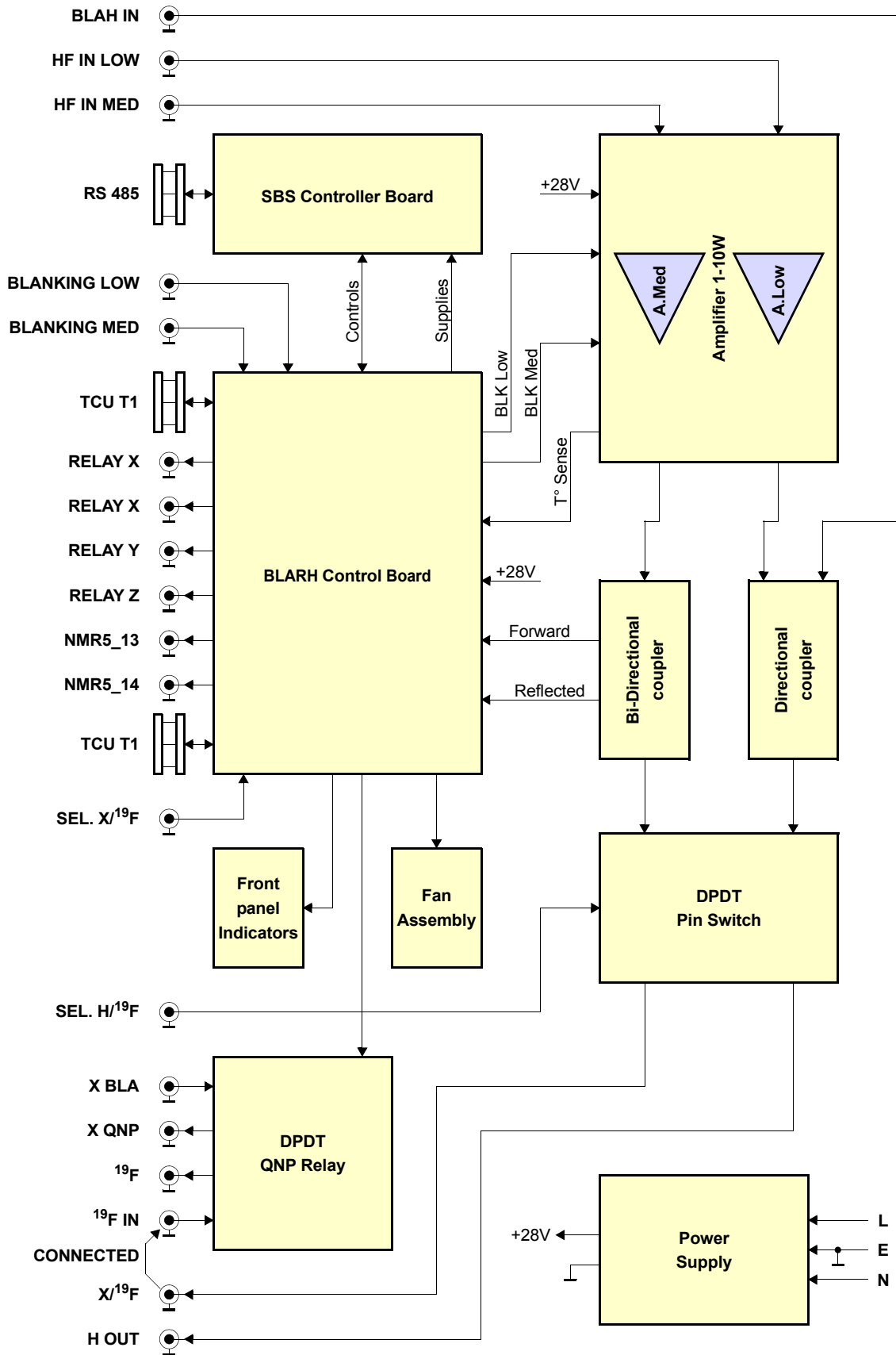

Technical description

5

System overview

5.1

The BLARH1-10 Amplifier 200-600MHz provide a RF Output Power of 10W and 1W in the 200-600MHz frequency range, for 1H or 19F channels on two distinct ways which can be routed via a DPDT switch to reach the front panel through 1H and 19F Outputs.

A special 1H input "BLAH IN" allows the mixing of an 100W external 1H ,like a BLAH100, with the 1W channel with the typical 50dB ratio.

With the "QNP Relay Option", the user can route via an internal DPDT relay an external X source, issued from a BALX300 for example, and the 10W channel to supply the both 19F and X QNP Outputs.

The linear module BLMH1-10 includes two class AB power amplifiers.

The 10W channel is connected to the DPDT Pin Switch via a Bi-Directional coupler. Reflect and Forward Information is conveyed to the user via the RS485 Bus on front panel.

The 1W channel is injected in the 100W 1H channel circuitry via a 50dB injection coupler.

The entire system is tied together by a Digital Signal Interface send power level monitor and malfunctions (supply and temperature faults) to the spectrometer.

Circuits such, Fan Fault or RF Module Overheating, Supply Status and 10W MED Channel activity are displayed on the Front Panel.

Theory of operation

5.2

RF Amplifier path


5.2.1

The BLARH1-10 Amplifier 200-600MHz (P/N: W1345023) consists of two Class AB power amplifier.

- A nominal input power level of +4dBm on RF «MED HF IN» produce a nominal output power of 10W peak or continuous power on MED channel.
- A nominal input power level of +4dBm on RF «LOW HF IN» produce a nominal output power of 1W peak or continuous power on the injection coupler or roughly 0dBm (-50dB / 100W) on the 1H / 19F high power section.

Technical description

Figure 5.1. Amplifier module block diagram

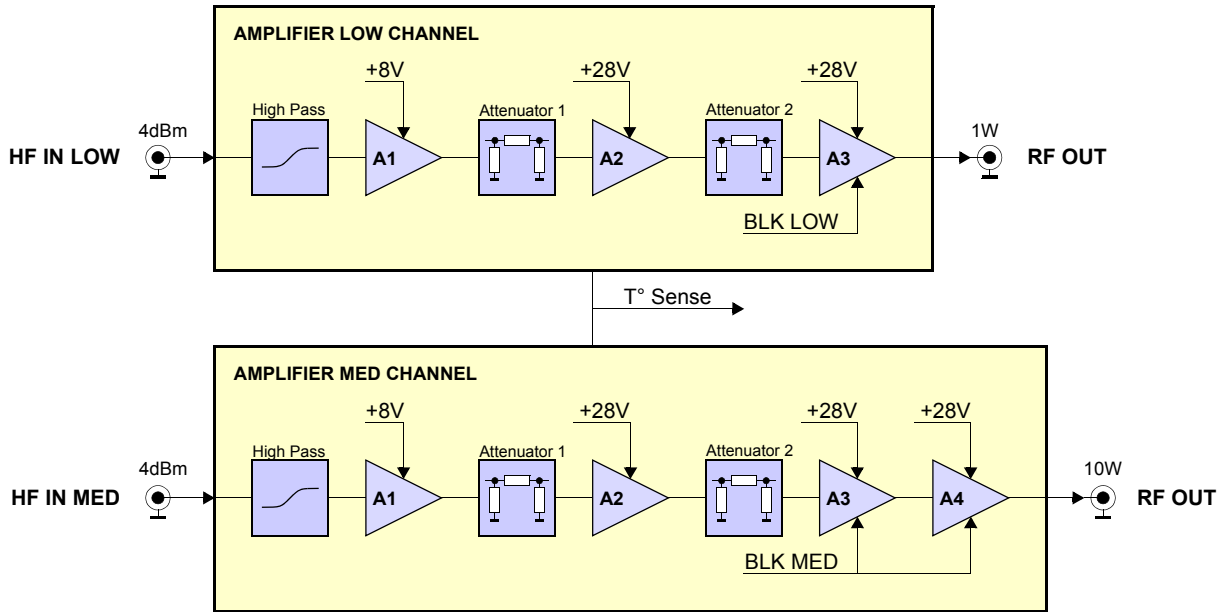

RF Outputs router without QNP Relay

5.2.2

Figure 5.2. DPDT Pin Switch routes


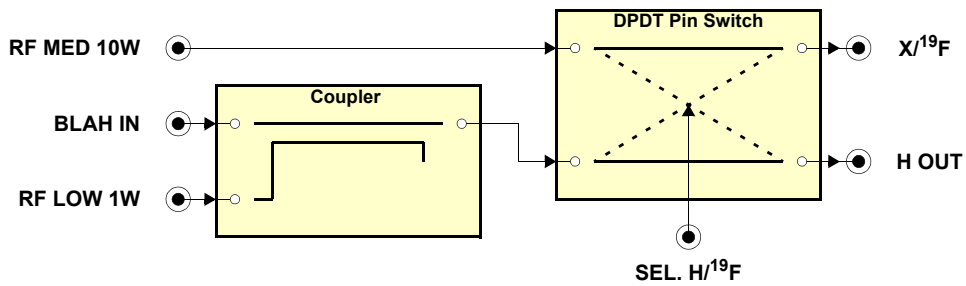


Table 5.1. DPDT Pin Switch True table

H/19F ^a	DPDT Pin Switch
0	Crossed
1	Through

a. TTL Level in BNC Connector.

Figure 5.3. View of ^{19}F Output

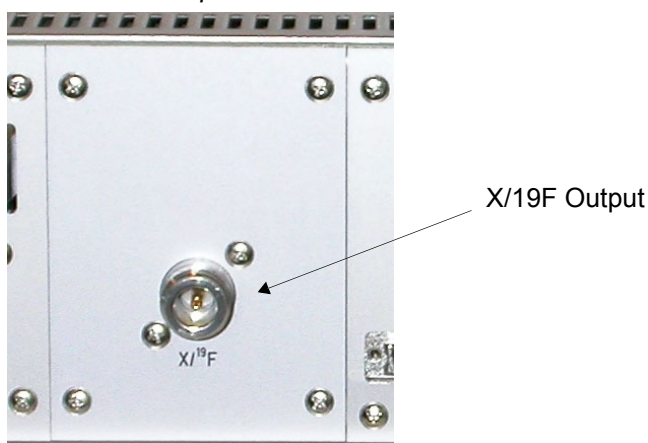

RF Outputs router with QNP Relay

5.2.3

Figure 5.4. DPDT Pin Switch and QNP Relay routes

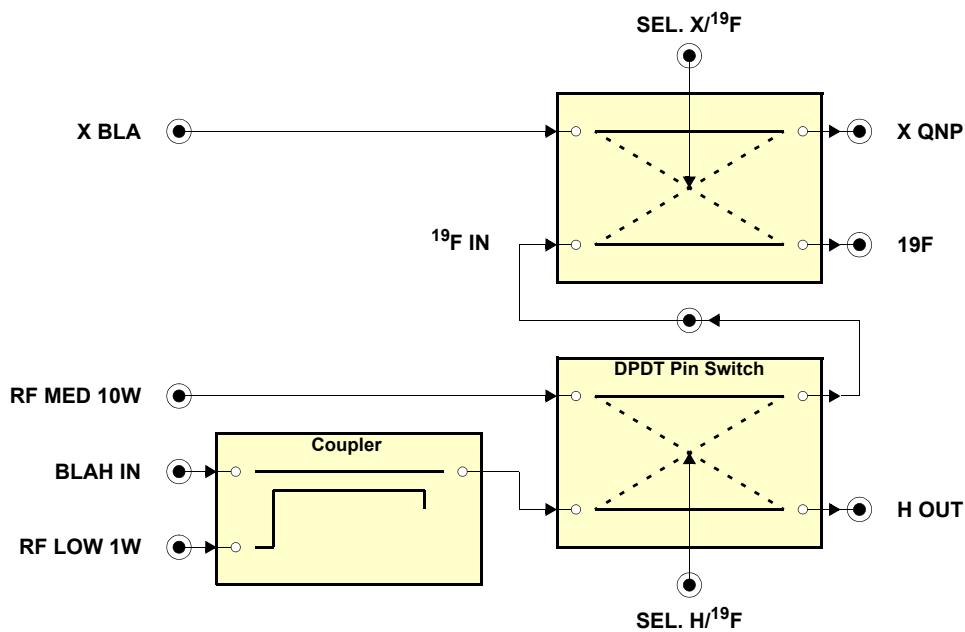


Table 5.2. DPDT Pin Switch

H/ ^{19}F ^a	DPDT Pin Switch
0	Crossed
1	Through

a. TTL Level in BNC Connector.

Table 5.3. DPDT QNP Relay True table

X/ ¹⁹ F ^a	DPDT QNP Relay
0	Crossed
1	Through

a. TTL Level in BNC Connector.

Figure 5.5. View of QNP Relay Accessory before installation



Control Board

5.2.4

The Control Board consists of circuitry to monitor the 10W output (MED Channel) RF power from the Bi-Directional coupler and to condition the input blanking (gating) signals and deliver there to the RF parts.

The Control Board also monitors the RF Path heatsink temperature and fans to protect against thermal overstress.

Information from supplies also being analyzed by the Control Board.

This Information are displayed on the front panel.

Status Board

5.2.5

The Status Board on the front panel of the amplifier displays overstress functions, supply status and "MED channel" (10W) RF activity.

Specifications

6

General specifications

6.1

Table 6.1. HMED and HLOW channel RF specifications

RF Specifications	Channel HMED	Channel HLOW
Frequency range	180 to 600MHz (3H on request)	180 to 600MHz (3H on request)
Linear Gain	38dB \pm 1dB	-2dB \pm 1dB
Gain Flatness	\pm 1dB max.	\pm 1,5dB max.
Minimum Pulsed Output Power	8W min.	-50dB of HHIGH Linear Power Region
CW Output Power	No limitation	No limitation
Linear Output Power	5W min. @ 1dBm compression	Full linear
Amplifier Biasing	Class AB Operation	Class A Operation
Blanking Delay	< 1 μ s typ.	< 1 μ s typ.
RF Rise Time	< 100ns	< 100ns
RF Fall Time	< 50ns	< 50ns
DC Ringing	200mV typ. (due to blanking signal)	100mV typ. (due to blanking signal)
Input Noise Figure	7dB max.	7dB max.
Output Noise Power (Unblanked)	-129dBm @ 1Hz	-169dBm @ 1Hz
Output Noise Power (Blanked)	-174dBm @ 1Hz (Thermal Noise)	-174dBm @ 1Hz (Thermal Noise)
IN/OUT Impedance	50 Ω	50 Ω
Input V.S.W.R.	1,5 max.	1,5 max.
Output Harmonics	30dBc min. (full range)	40dBc min. (full range)
Pulse Width (int. limitation)	No limitation	No limitation
Duty Cycle (int. limitation)	No limitation	No limitation
Amplitude Droop	< 4% @ 10W for 500ms Pulse Width	< 1% (full power, full range)

Table 6.2. External H150 channel RF specifications

RF Specifications	Channel H150
Frequency range	188 to 400MHz / 470 to 500MHz (3H on request)
Linear Gain	48dB ±1dB
Gain Flatness	±2dB max.
Minimum Pulsed Output Power	150W typ. (@ nominal input +4dBm)
CW Output Power	25W max. (internal limitation)
Linear Output Power	120W min. @ 1dBm compression
Amplifier Biasing	Class AB Operation
Blanking Delay	< 1µs typ.
RF Rise Time	< 70ns
RF Fall Time	< 50ns
Input Noise Figure	6dB max.
Output Noise Power (Unblanked)	-118dBm @ 1Hz
Output Noise Power (Blanked)	-174dBm @ 1Hz (Thermal Noise)
IN/OUT Impedance	50Ω
Input V.S.W.R.	1,5 max.
Output Harmonics	-30dBc min. @ 150W (full range)
Pulse Width (int. limitation)	250ms @ 200W (up to CW @ 50W)
Duty Cycle (int. limitation)	5% @ 200W (up to 100% @ 50W)
Amplitude Droop	< 6% @ 200W for 100ms Pulse Width
Dimensions	19" rack cabinet, 2U height, 460mm deph
Weight	12kg

RF signal coming from BLAH1000 (High Resolution output). Signal controlled on High Power "H" Amplifier (coupled under 50dB with HLOW channel). Internal insertion loss : 2dB.

Table 6.3. Common characteristics

RF Input	SMA (F)
RF Output	N (F)
Blanking Pulse Connector	BNC (F)
RF Switch Control Signal	BNC (F)

Figures

1	General information	5
Figure 1.1.	BLARH1-10 Amplifier 200-600MHz without XQNP accessory	5
2	Safety	7
Figure 2.1.	Identifying plate	7
Figure 2.2.	Manufacturer's nameplate	8
Figure 2.3.	Electrical hazard symbol	9
Figure 2.4.	General hazard symbol	9
3	Installation	11
4	Operation	13
Figure 4.1.	BLARH1-10 Front Panel Design without QNP option	13
Figure 4.2.	BLARH1-10 Front Panel Design with QNP option	13
Figure 4.3.	Amplifier block diagram without QNP option	18
Figure 4.4.	Amplifier block diagram with QNP option	19
5	Technical description	21
Figure 5.1.	Amplifier module block diagram	22
Figure 5.2.	DPDT Pin Switch routes	22
Figure 5.3.	View of ¹⁹ F Output	23
Figure 5.4.	DPDT Pin Switch and QNP Relay routes	23
Figure 5.5.	View of QNP Relay Accessory before installation	24
6	Specifications	25

Figures

Tables

1	General information	5
2	Safety	7
3	Installation	11
4	Operation	13
Table 4.1.	LED's indicators assignment	14
Table 4.2.	RF Connectors assignment without QNP Relay Option ...	14
Table 4.3.	RF Connectors assignment with QNP Relay Option	14
Table 4.4.	Blanking Connectors assignment	15
Table 4.5.	RF Output Selection	15
Table 4.6.	DPDT Pin Switch routing true table	15
Table 4.7.	DPDT QNP Relay routing true table	15
Table 4.8.	TCU T1 Connectors assignment	16
Table 4.9.	Correspondance between TCU T1 and SMB Connectors	17
5	Technical description	21
Table 5.1.	DPDT Pin Switch True table	22
Table 5.2.	DPDT Pin Switch	23
Table 5.3.	DPDT QNP Relay True table	24
6	Specifications	25
Table 6.1.	HMED and HLOW channel RF specifications	25
Table 6.2.	External H150 channel RF specifications	26
Table 6.3.	Common characteristics	26

